



**APEX BIO**  
Achieve Perfection, Explore the Unknown

# High Performance Liquid Chromatography


## 1. Analysis Information

Product Name:	Gap 26
Operator:	David
Injection Date:	11/24/2017 3:47:59 PM
Batch No.:	2

## 2. HPLC Condition

Column:	Athena C18, 3 $\mu$ m, 2.1mm $\times$ 100 mm																			
Solvent A:	0.1% H3PO4 in 100% Acetonitrile																			
Solvent B:	0.1% H3PO4 in 100% Water																			
Gradient:	<table border="1"> <thead> <tr> <th>Time</th> <th>A</th> <th>B</th> </tr> </thead> <tbody> <tr> <td>0.0min</td> <td>0%</td> <td>100%</td> </tr> <tr> <td>2.0min</td> <td>0%</td> <td>100%</td> </tr> <tr> <td>7.0min</td> <td>95%</td> <td>5%</td> </tr> <tr> <td>10.0 min</td> <td>95%</td> <td>5%</td> </tr> <tr> <td>10.0min</td> <td>Stop</td> <td>Stop</td> </tr> </tbody> </table>		Time	A	B	0.0min	0%	100%	2.0min	0%	100%	7.0min	95%	5%	10.0 min	95%	5%	10.0min	Stop	Stop
Time	A	B																		
0.0min	0%	100%																		
2.0min	0%	100%																		
7.0min	95%	5%																		
10.0 min	95%	5%																		
10.0min	Stop	Stop																		
Flow rate:	0.45ml/min																			
Wavelength:	220nm																			
Volume:	10 $\mu$ l																			

## 3. Result

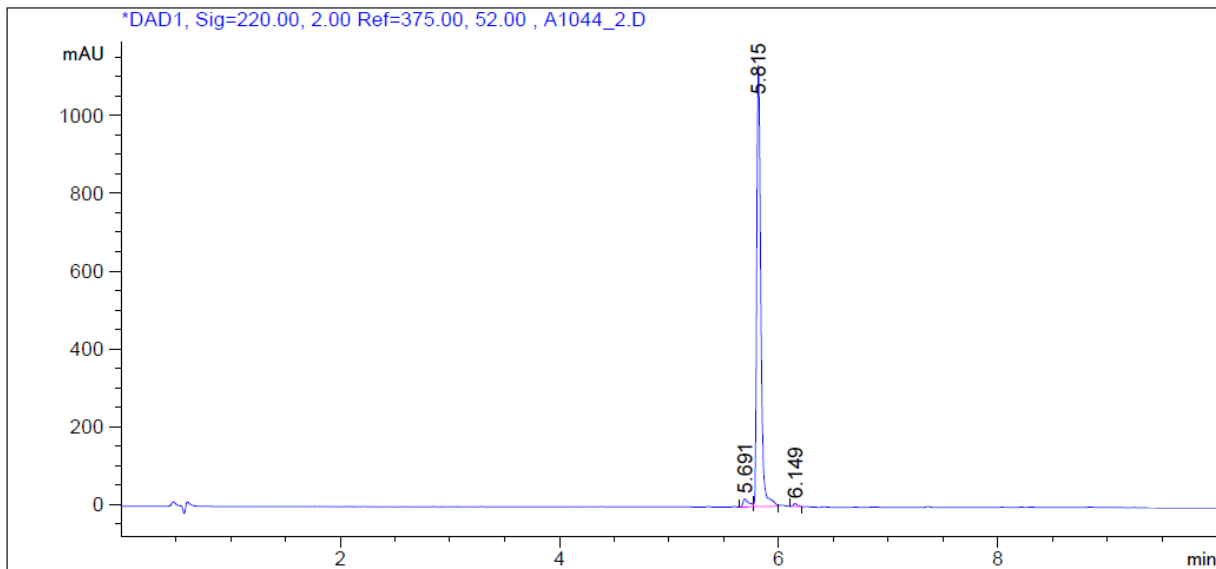

Rank	RetTime(min)	Area(mAU*s)	Area %
1	5.691	45.96255	1.5434
2	5.815	2919.53223	98.0361
3	6.149	12.52207	0.4205