

High Performance Liquid Chromatography

Analysis Information

Product Name: MEK162 (ARRY-162, ARRY-438162)
Operator: David
Injection Date: 7/12/2019 9:47:38 PM
Batch No.: 2

HPLC Condition

Column : ZORBAX Eclipse Pluse C18 3.5 μ m 4.6mmx100mm
Solvent A: H2O+1‰ TFA
Solvent B: ACN+1‰ TFA
Gradient:

Time	A	B
0.0min	80%	20%
2.0min	80%	20%
7.0min	5%	95%
10.0min	0%	100%
12.0min	0%	100%
12.1min	80%	20%
15.0min	Stop	Stop


Flow rate: 1ml/min

Wavelength: 250nm

Volume:

5 μ l

Result

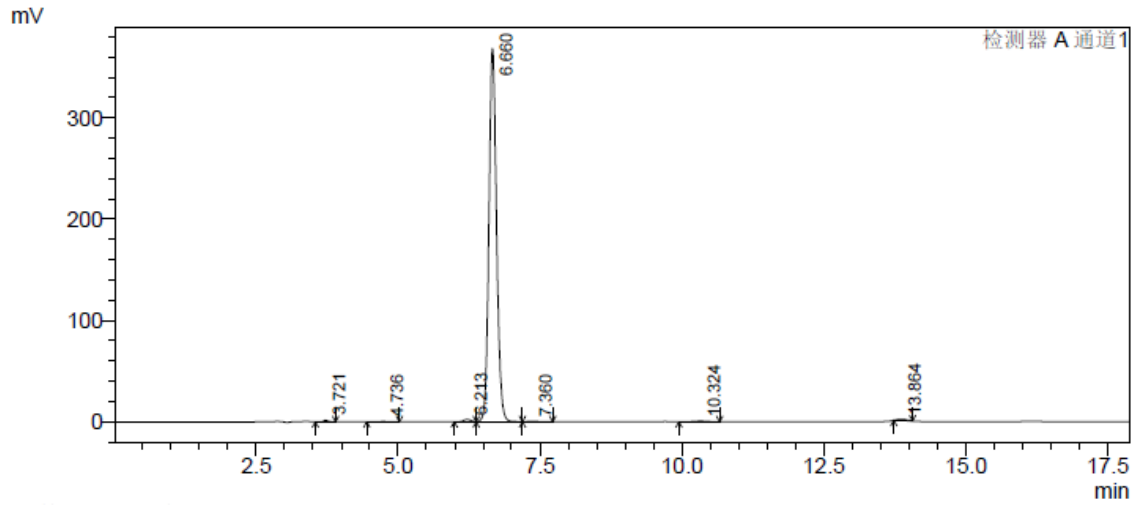

Rank	RetTime(min)	Area(mAU*s)	Area %
1	3.721	7976	0.370
2	4.736	3462	0.096
3	6.213	21494	0.605
4	6.660	3441793	98.178
5	7.360	6330	0.150
6	10.324	11902	0.217
7	13.864	16893	0.384

Caution

FOR RESEARCH PURPOSES ONLY.

NOT FOR HUMAN, VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

APEXIO Technology

www.apexbt.com

7505 Fannin street, Suite 410, Houston, TX 77054.

Tel: +1-832-696-8203 | Fax: +1-832-641-3177 | Email: info@apexbt.com