

High Performance Liquid Chromatography

Analysis Information

Product Name: NVP-BSK805 2HCl
Operator: David
Injection Date: 4/22/2020 9:09:37 AM
Batch No.: 2

HPLC Condition


Column : ZORBAX Eclipse Pluse C18 3.5 μ m 4.6mmx100mm
Solvent A: H2O+0.1%TFA
Solvent B: ACN+0.1%TFA

Gradient:

| Time | A | B |
|----------|------|------|
| 0.0min | 80% | 20% |
| 2.0min | 80% | 20% |
| 7.0min | 5% | 95% |
| 10.00min | 0% | 100% |
| 12.00min | 0% | 100% |
| 15.00min | Stop | Stop |

Flow rate: 1mL/min
Wavelength: 280nm
Volume: 5 μ l

Result

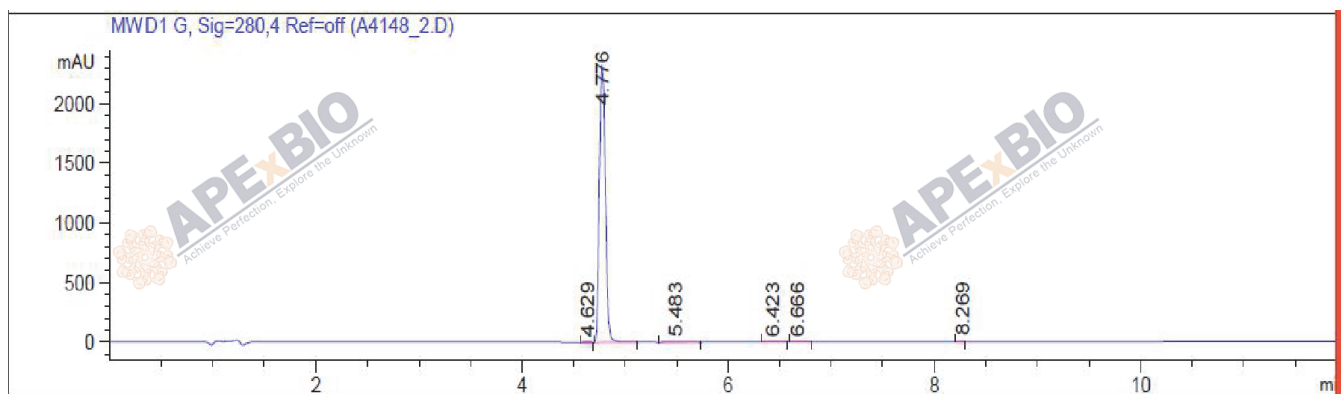

| Rank | RetTime(min) | Area(mAU*s) | Area % |
|------|--------------|-------------|---------|
| 1 | 4.629 | 20.34808 | 0.2144 |
| 2 | 4.776 | 9418.91406 | 99.2341 |
| 3 | 5.483 | 16.20782 | 0.1708 |
| 4 | 6.423 | 16.21229 | 0.1708 |
| 5 | 6.666 | 10.98119 | 0.1157 |
| 6 | 8.269 | 8.94709 | 0.0943 |

Caution

FOR RESEARCH PURPOSES ONLY.

NOT FOR HUMAN, VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

APEX BIO Technology

www.apexbt.com

7505 Fannin street, Suite 410, Houston, TX 77054.

Tel: +1-832-696-8203 | Fax: +1-832-641-3177 | Email: info@apexbt.com