

# High Performance Liquid Chromatography

## Analysis Information

**Product Name:** SAR407899  
**Operator:** David  
**Injection Date:** 12/14/2023 11:29:05 AM  
**Batch No.:** 2


## HPLC Condition

**Column :** ZORBAX Eclipse Pluse C18 3.5 $\mu$ m 4.6mmx100mm  
**Solvent A:** H2O+0.1%TFA  
**Solvent B:** ACN+0.1%TFA  
**Gradient:**

| Time    | A    | B    |
|---------|------|------|
| 0.0min  | 80%  | 20%  |
| 2.0min  | 80%  | 20%  |
| 7.0min  | 5%   | 95%  |
| 10.0min | 0%   | 100% |
| 12.0min | 0%   | 100% |
| 15.0min | Stop | Stop |

**Flow rate:** 1ml/min  
**Wavelength:** 280nm  
**Volume:** 2 $\mu$ l

# Result

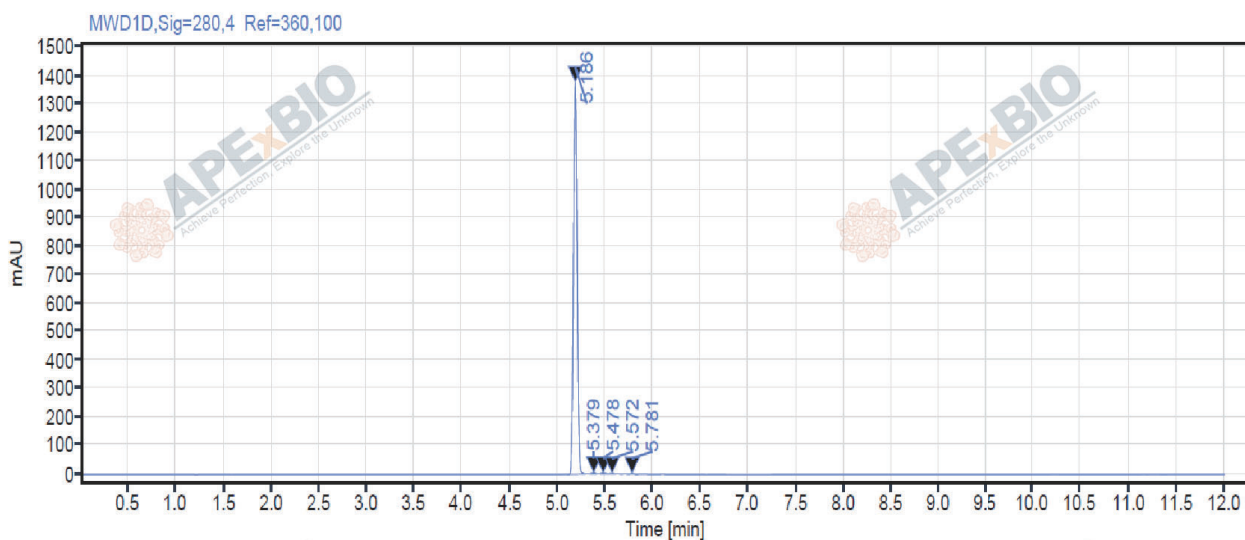

| Rank | RetTime(min) | Area(mAU*s) | Area % |
|------|--------------|-------------|--------|
| 1    | 5.186        | 3776.05     | 99.10  |
| 2    | 5.379        | 10.19       | 0.27   |
| 3    | 5.478        | 10.26       | 0.27   |
| 4    | 5.572        | 6.99        | 0.18   |
| 5    | 5.781        | 6.84        | 0.18   |

## Caution

**FOR RESEARCH PURPOSES ONLY.**

**NOT FOR HUMAN, VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.**

**APEX BIO Technology**

[www.apexbt.com](http://www.apexbt.com)

7505 Fannin Street, Suite 410, Houston, TX 77054.

Tel: +1-832-696-8203 | Fax: +1-832-641-3177 | Email: [info@apexbt.com](mailto:info@apexbt.com)

