

# High Performance Liquid Chromatography

## Analysis Information

**Product Name:** MDV3100 (Enzalutamide)

**Operator:** David

**Injection Date:** 11/29/2019 6:43:22 PM

**Batch No.:** 4

## HPLC Condition

**Column :** ZORBAX Eclipse Pluse C18 3.5 $\mu$ m 4.6mmx100mm

**Solvent A:** H2O+1‰ TFA

**Solvent B:** ACN+1‰ TFA

**Gradient:**


| Time    | A    | B    |
|---------|------|------|
| 0.0min  | 80%  | 20%  |
| 2.0min  | 80%  | 20%  |
| 7.0min  | 5%   | 95%  |
| 10.0min | 0%   | 100% |
| 12.0min | 0%   | 100% |
| 12.1min | 80%  | 20%  |
| 15.0min | Stop | Stop |

**Flow rate:** 1ml/min

**Wavelength:** 280nm

**Volume:** 5 $\mu$ l

Result

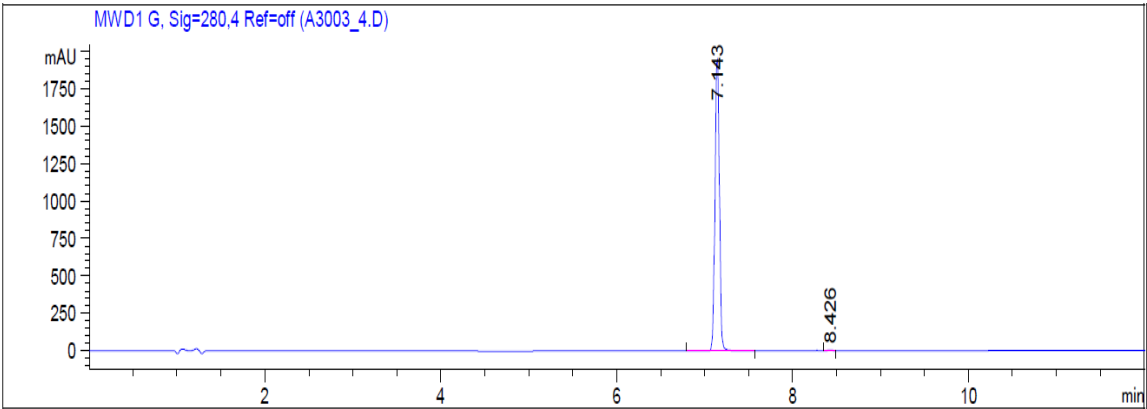

| Rank | RetTime(min) | Area(mAU*s) | Area %  |
|------|--------------|-------------|---------|
| 1    | 7.143        | 6812.29785  | 99.7845 |
| 2    | 8.426        | 14.71101    | 0.2155  |

Caution

FOR RESEARCH PURPOSES ONLY.

NOT FOR HUMAN, VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

APExBIO Technology

www.apexbt.com

7505 Fannin street, Suite 410, Houston, TX 77054.

Tel: +1-832-696-8203 | Fax: +1-832-641-3177 | Email: info@apexbt.com