

# High Performance Liquid Chromatography

## Analysis Information

**Product Name:** TBB

**Operator:** David

**Injection Date:** 4/22/2020 12:19:59 AM

**Batch No.:** 2

## HPLC Condition

**Column :** ZORBAX Eclipse Pluse C18 3.5 $\mu$ m 4.6mmx100mm

**Solvent A:** H2O+0.1%TFA

**Solvent B:** ACN+0.1%TFA

**Gradient:**


| Time     | A    | B    |
|----------|------|------|
| 0.0min   | 80%  | 20%  |
| 2.0min   | 80%  | 20%  |
| 7.0min   | 5%   | 95%  |
| 10.00min | 0%   | 100% |
| 12.00min | 0%   | 100% |
| 15.00min | Stop | Stop |

**Flow rate:** 1mL/min

**Wavelength:** 270nm

**Volume:** 5 $\mu$ l

## Result

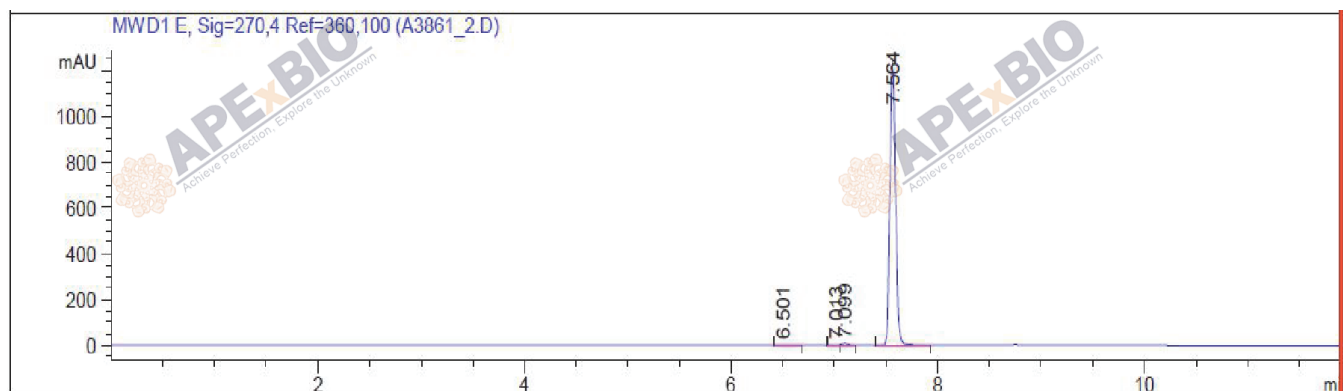

| Rank | RetTime(min) | Area(mAU*s) | Area %  |
|------|--------------|-------------|---------|
| 1    | 6.501        | 4.06857     | 0.0847  |
| 2    | 7.013        | 3.40829     | 0.0709  |
| 3    | 7.099        | 29.19120    | 0.6077  |
| 4    | 7.564        | 4767.28076  | 99.2367 |

## Caution

**FOR RESEARCH PURPOSES ONLY.**

**NOT FOR HUMAN, VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.**

**APEX BIO Technology**

**[www.apexbt.com](http://www.apexbt.com)**

7505 Fannin street, Suite 410, Houston, TX 77054.

Tel: +1-832-696-8203 | Fax: +1-832-641-3177 | Email: [info@apexbt.com](mailto:info@apexbt.com)