

# High Performance Liquid Chromatography

## Analysis Information

**Product Name:** Paliperidone  
**Operator:** David  
**Injection Date:** 3/26/2019 10:00:06 PM  
**Batch No.:** 1


## HPLC Condition

**Column :** ZORBAX Eclipse Pluse C18 3.5 $\mu$ m 4.6mmx100mm  
**Solvent A:** H2O+1‰ TFA  
**Solvent B:** ACN+1‰ TFA  
**Gradient:**

| Time    | A    | B    |
|---------|------|------|
| 0.0min  | 95%  | 5%   |
| 2.0min  | 95%  | 5%   |
| 7.0min  | 5%   | 95%  |
| 10.0min | 0%   | 100% |
| 12.0min | 0%   | 100% |
| 12.1min | 95%  | 5%   |
| 15.0min | Stop | Stop |

**Flow rate:** 1ml/min  
**Wavelength:** 280nm  
**Volume:** 5 $\mu$ l

Result

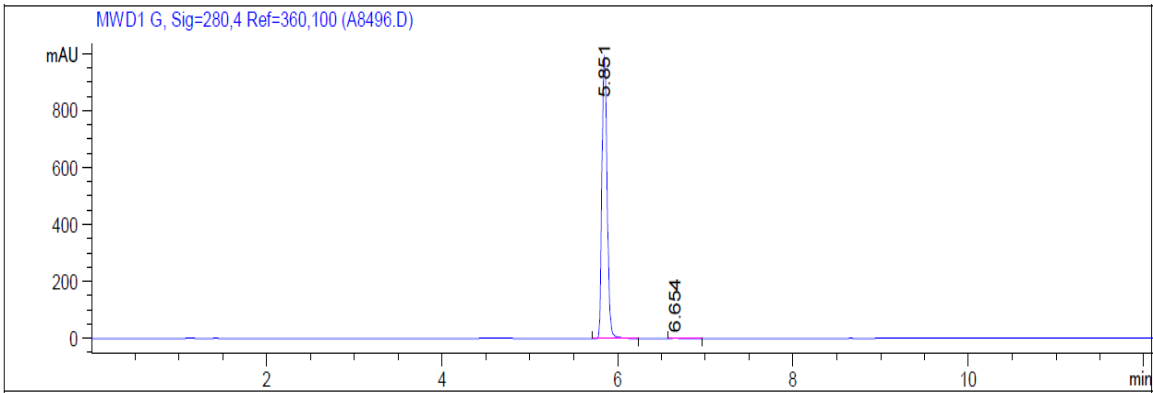

| Rank | RetTime(min) | Area(mAU*s) | Area %  |
|------|--------------|-------------|---------|
| 1    | 5.851        | 4124.93555  | 99.8161 |
| 2    | 6.654        | 7.59985     | 0.1839  |

Caution

FOR RESEARCH PURPOSES ONLY.

NOT FOR HUMAN, VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

APExBIO Technology

www.apexbt.com

7505 Fannin street, Suite 410, Houston, TX 77054.

Tel: +1-832-696-8203 | Fax: +1-832-641-3177 | Email: info@apexbt.com