

High Performance Liquid Chromatography

Analysis Information

Product Name: Safe DNA Gel Stain
Operator: David
Injection Date: 2021/12/20 14:20:06
Batch No.: 8


HPLC Condition

Column : ZORBAX Eclipse Pluse C18 3.5 μ m 4.6mmx100mm
Solvent A: H2O+0.1%TFA
Solvent B: ACN+0.1%TFA
Gradient:

Time	A	B
0.0min	80%	20%
2.0min	80%	20%
7.0min	5%	95%
10.0min	0%	100%
12.0min	0%	100%
15.0min	Stop	Stop

Flow rate: 1ml/min
Wavelength: 220nm
Volume: 5 μ l

Result

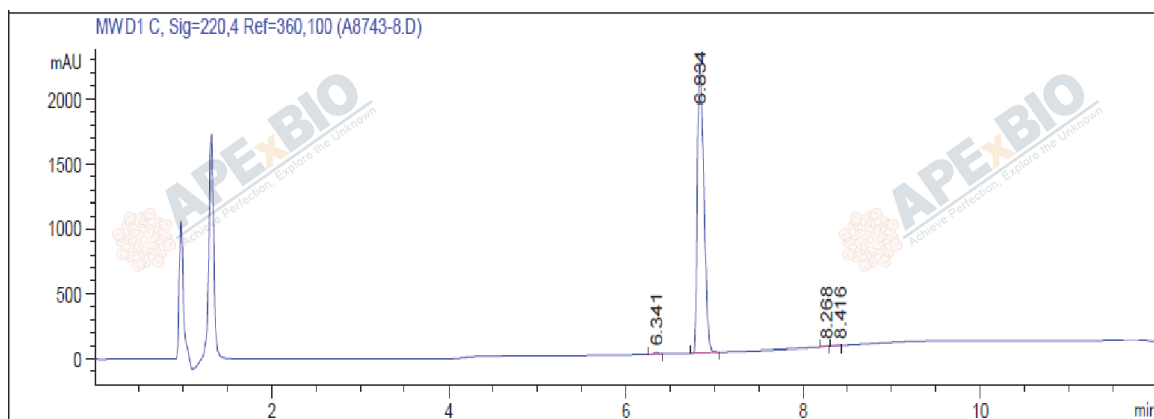

Rank	RetTime(min)	Area(mAU*s)	Area %
1	6.341	43.79019	0.3717
2	6.834	1.17246e4	99.5119
3	8.268	5.42157	0.0460
4	8.416	8.29655	0.0704

Caution

FOR RESEARCH PURPOSES ONLY.

NOT FOR HUMAN, VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

APEX BIO Technology

www.apexbt.com

7505 Fannin street, Suite 410, Houston, TX 77054.

Tel: +1-832-696-8203 | Fax: +1-832-641-3177 | Email: info@apexbt.com