

High Performance Liquid Chromatography

Analysis Information

Product Name: Pyridostatin TFA

Operator: David

Injection Date: 07/09/2024 13:06:36

Batch No.: 2

HPLC Condition

Column : ZORBAX Eclipse Pluse C18 3.5 μ m 4.6mmx100mm

Solvent A: H2O+0.1%TFA

Solvent B: ACN+0.1%TFA

Gradient:


Time	A	B
0.0min	95%	5%
2.0min	95%	5%
16.0min	5%	95%
17.0min	0%	100%
20.0min	0%	100%
23.0min	Stop	Stop

Flow rate: 1ml/min

Wavelength: 254nm

Volume: 5 μ l

Result

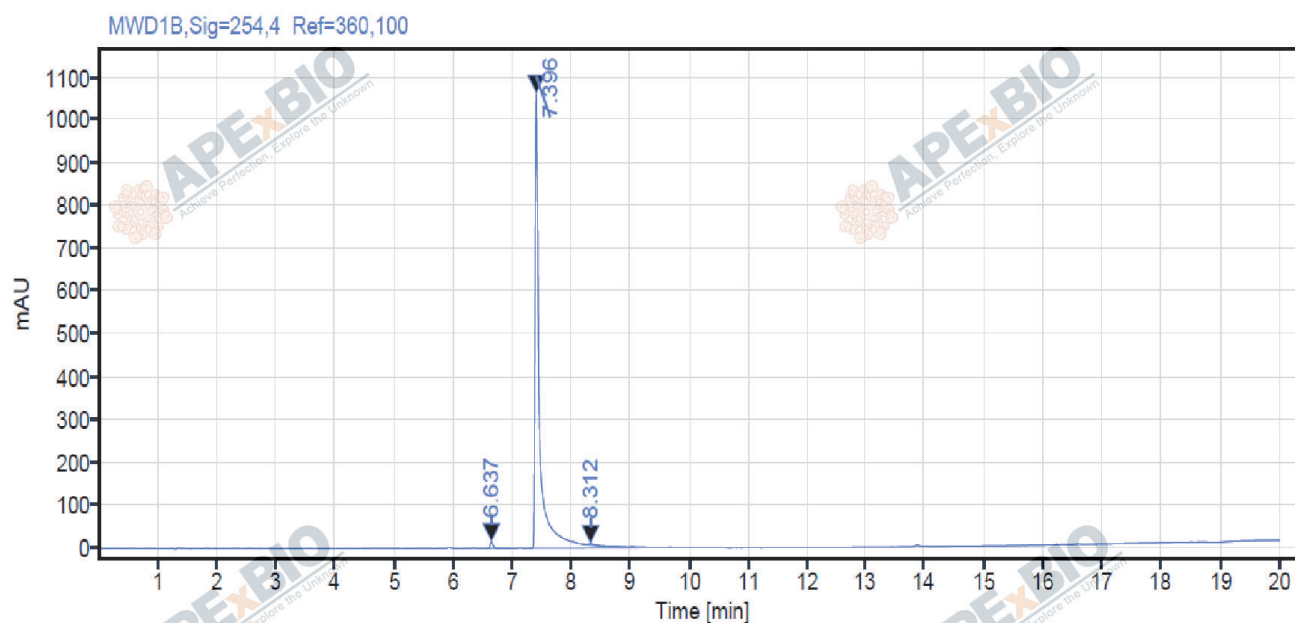

Rank	RetTime(min)	Area(mAU*s)	Area %
1	6.637	54.12	0.92
2	7.396	5782.21	98.22
3	8.312	50.82	0.86

Caution

FOR RESEARCH PURPOSES ONLY.

NOT FOR HUMAN, VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

APExBIO Technology

www.apexbt.com

7505 Fannin street, Suite 410, Houston, TX 77054.

Tel: +1-832-696-8203 | Fax: +1-832-641-3177 | Email: info@apexbt.com