

High Performance Liquid Chromatography

Analysis Information

Product Name: HSP990 (NVP-HSP990)

Operator: David

Injection Date: 6/10/2020 7:42:56 PM

Batch No.: 2

HPLC Condition

Column : ZORBAX Eclipse Pluse C18 3.5 μ m 4.6mmx100mm

Solvent A: H2O+0.1%TFA

Solvent B: ACN+0.1%TFA

Gradient:


Time	A	B
0.0min	80%	20%
2.0min	80%	20%
7.0min	5%	95%
10.0min	0%	100%
12.0min	0%	100%
15.0min	Stop	Stop

Flow rate: 1ml/min

Wavelength: 280nm

Volume: 5 μ l

Result

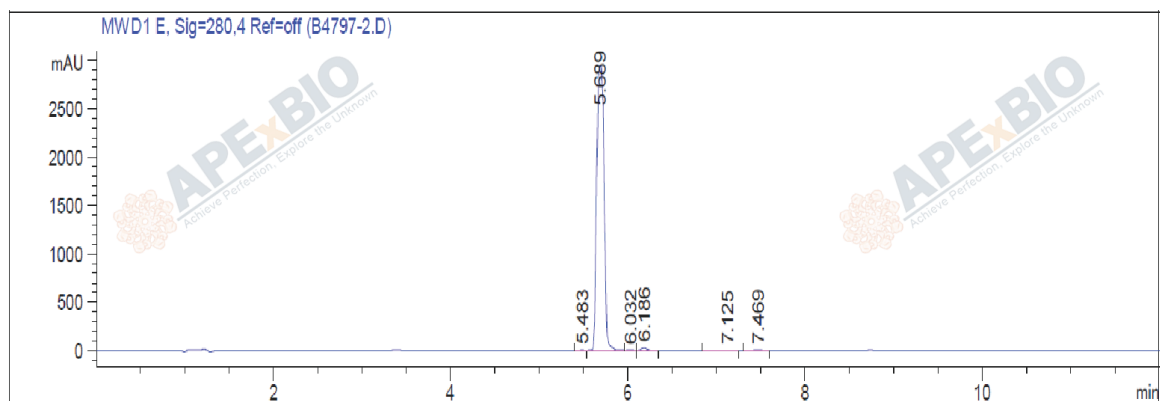

Rank	RetTime(min)	Area(mAU*s)	Area %
1	5.483	10.90118	0.0650
2	5.689	1.66016e4	99.0187
3	6.032	11.81781	0.0705
4	6.186	111.56044	0.6654
5	7.125	6.83562	0.0408
6	7.469	23.41980	0.1397

Caution

FOR RESEARCH PURPOSES ONLY.

NOT FOR HUMAN, VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

APExBIO Technology

www.apexbt.com

7505 Fannin street, Suite 410, Houston, TX 77054.

Tel: +1-832-696-8203 | Fax: +1-832-641-3177 | Email: info@apexbt.com