

High Performance Liquid Chromatography

Analysis Information

Product Name: (±)-Nutlin-3

Operator: David

Injection Date: 03/20/2024 14:30:24

Batch No.: 2

HPLC Condition

Column : ZORBAX Eclipse Pluse C18 3.5 μ m 4.6mmx100mm

Solvent A: H2O+0.1%TFA

Solvent B: ACN+0.1%TFA

Gradient:


Time	A	B
0.0min	80%	20%
2.0min	80%	20%
7.0min	5%	95%
10.0min	0%	100%
12.0min	0%	100%
15.0min	Stop	Stop

Flow rate: 1ml/min

Wavelength: 280nm

Volume: 2 μ l

Result

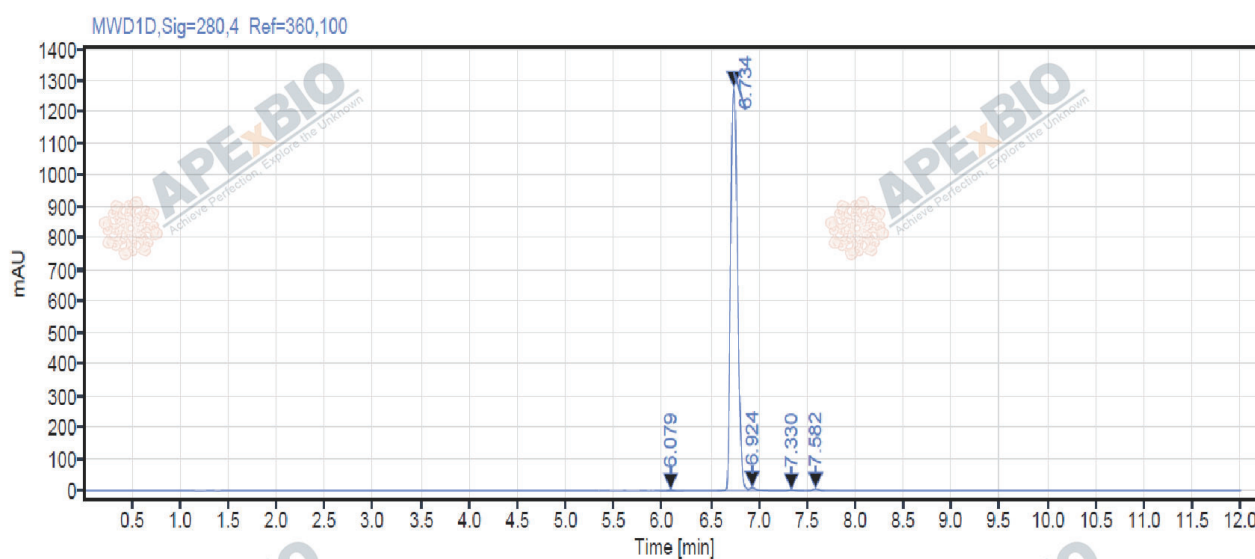

Rank	RetTime(min)	Area(mAU*s)	Area %
1	6.079	12.48	0.19
2	6.734	6441.01	98.45
3	6.924	53.07	0.81
4	7.330	11.69	0.18
5	7.582	24.32	0.37

Caution

FOR RESEARCH PURPOSES ONLY.

NOT FOR HUMAN, VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

APExBIO Technology

www.apexbt.com

7505 Fannin street, Suite 410, Houston, TX 77054.

