

High Performance Liquid Chromatography

Analysis Information

Product Name: proTAME

Operator: David

Injection Date: 02/14/2025 13:19:20

Batch No.: 2

HPLC Condition

Column : ZORBAX Eclipse Pluse C18 3.5 μ m 4.6mmx100mm

Solvent A: H2O+0.1%TFA

Solvent B: ACN+0.1%TFA

Gradient:


Time	A	B
0.0min	80%	20%
2.0min	80%	20%
7.0min	5%	95%
10.0min	0%	100%
12.0min	0%	100%
15.0min	Stop	Stop

Flow rate: 1ml/min

Wavelength: 220nm

Volume: 0.3 μ l

Result

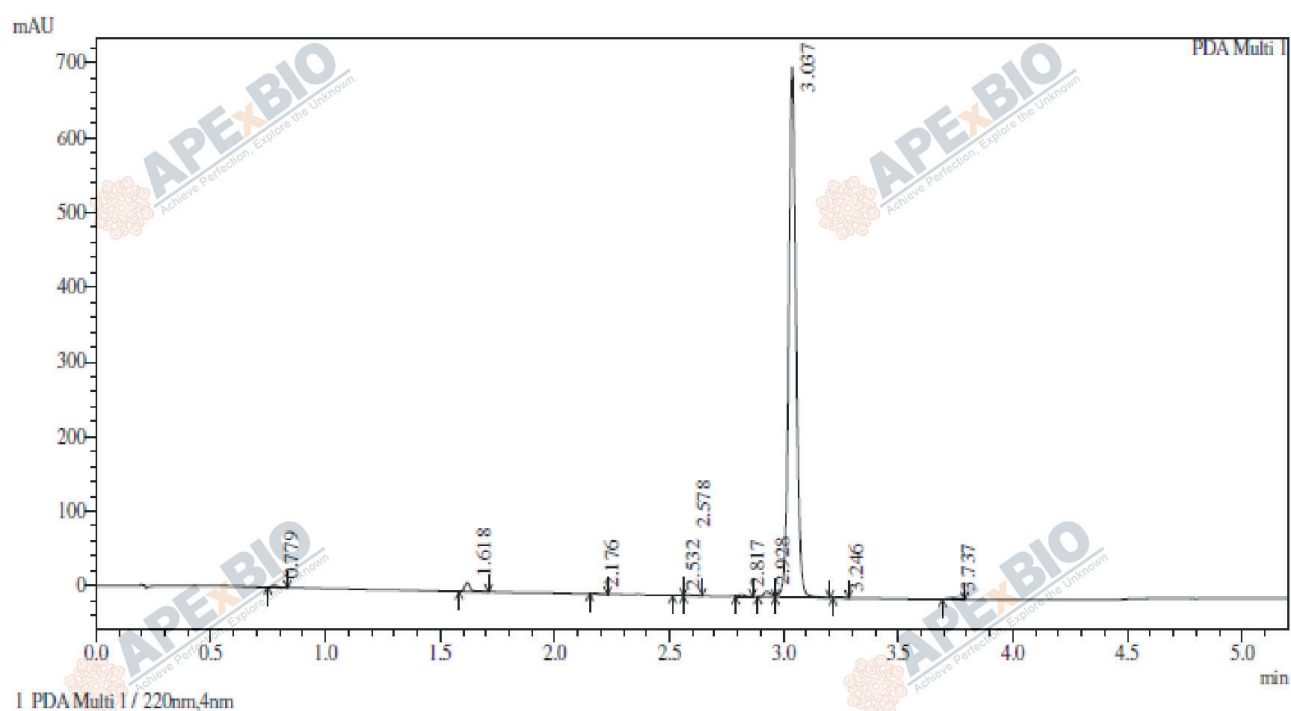

Rank	RetTime(min)	Area(mAU*s)	Area %
1	0.779	3333	0.204
2	1.618	19039	1.168
3	2.176	1546	0.095
4	2.532	1168	0.072
5	2.578	2829	0.174
6	2.817	4849	0.297
7	2.928	15197	0.932
8	3.037	1575324	96.635
9	3.246	1770	0.109
10	3.737	5126	0.314

Caution

FOR RESEARCH PURPOSES ONLY.

NOT FOR HUMAN, VETERINARY DIAGNOSTIC OR THERAPEUTIC USE.

APEx BIO Technology

www.apexbt.com

7505 Fannin street, Suite 410, Houston, TX 77054.

Tel: +1-832-696-8203 | Fax: +1-832-641-3177 | Email: info@apexbt.com

